



LEADERS HOLDING CO.™

Est. 1971

Scanning Temperature Touchless Thermometer

All Leaders Holding Co. facilities open to the public have been provided a touchless infrared thermometer to improve the safety of our team and guests. It is not legally required to scan team members or guests, but our teams now have the ability if the need arises.

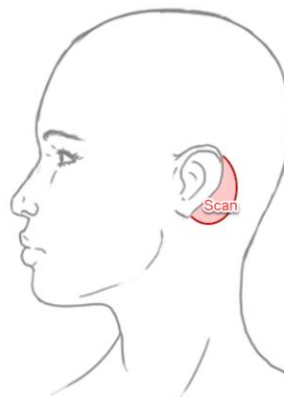
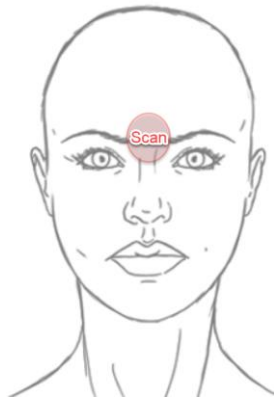
Temperature Readings

Please make sure your device is set to “Body.” See instructions on page 2 if you need help

	Fahrenheit	Celsius		Access to Facility
Safe	<99	<32	Normal temperature	Yes
Caution	99 - 100	37.4 – 38	Slightly elevated. Could be from activity, the sun, or just coming in from outside on a warm day	Measure Multiple Times, Store Teams Decision
Danger	101 <	38.1 <	Fever very likely	No



How to Scan

- ✓ Make sure eyeglasses and hair are out of the way
- ✓ If a person is perspiring, use the behind the ear measurement
- ✓ From one to two inches back, aim the device at the forehead between the eyebrows
- ✓ Take more than one reading for accuracy



Device Instructions

Leaders has multiple devices in circulation. Please look at the device instruction online for specific operating instructions.

Device	Model	Link to manual
	ZS-T1	<p>How to turn the power on / off:</p> <ol style="list-style-type: none">1. Press the M button to turn on the power.2. Press the M button again to turn off the power <p>Automatic power-off without operation within 60 seconds</p> <p>Measurement mode switch:</p> <p>Press "-" button to switch between surface temperature mode and body temperature mode</p> <p>Switch between Celsius and Fahrenheit:</p> <p>Press and hold the "-" button for 5 seconds, and the F1 symbol will appear on the screen.</p> <p>Then press the "-" button three times to display the F4 symbol on the screen.</p> <p>Press the "+" button. You should see the temperature units change to Fahrenheit.</p> <p>Turn off power to save settings</p>
	NUB838OH	<p>How to turn the power on / off:</p> <ol style="list-style-type: none">1. Squeeze the trigger <p>Automatic power-off without operation within 60 seconds</p> <p>Measurement mode switch:</p> <p>Press the "MODE" button to switch between surface temperature mode and body temperature mode</p> <p>Switch between Celsius and Fahrenheit:</p> <p>Press and hold the "MODE" button for 5 seconds</p> <p>Press Up Arrow for Celsius</p> <p>Press Down Arrow for Fahrenheit</p> <p>Press and hold "MODE" button for 5 seconds to save</p>